

# Introducing the PUNCHTac™ Anchor System with Fiber-Connect™ Technology

Smaller Anchor, Stronger Fixation, and Easier Suture-Tape Passage.


Available in 4.5mm, 4.75mm and 5.5mm




NEW

## What is Fiber-Connect™ Technology?

Fiber-Connect™ Technology transitions from tape to a suture configuration to ease with passage through the tissue being repaired and assist with knot tying.

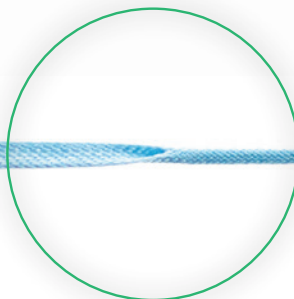

Fiber-Connect tape also provides a broader footprint at the site of repair as compared to a standard #2 suture.

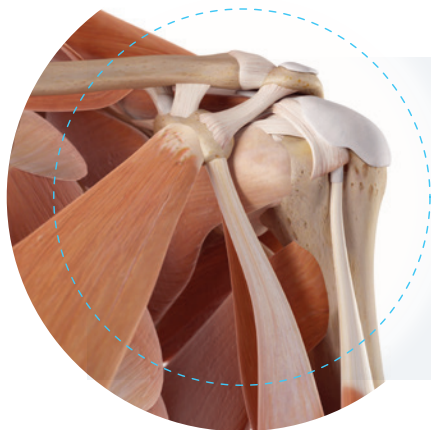

**The PUNCHTac™ 5.5mm** Fully Threaded Anchor introduces an advanced thread pitch design for the shoulder that's a game changer for medial row repair. It's about maximizing bone-to-anchor surface contact and distributing stress or torsional load evenly. All of this results in a Failure Force of 624.6 N<sup>1</sup>.

**DUNAMIS** MEDICAL™

Make the switch to better design

334.371.6280  
dunamismedical.com

Easier passage through soft-tissue



Effortless loading of knot pusher as compared to a standard tape




## More bone preservation without compromising fixation strength

### 4.5mm and 4.75mm Threaded PUNCHTac Anchor

The dual-lead thread design of the 4.5mm and 4.75mm anchors provide ease of insertion, enhanced fixation, and greater bone preservation at the time of implantation.



Average Load-to-Failure<sup>1</sup> - Soft Tissue Anchors

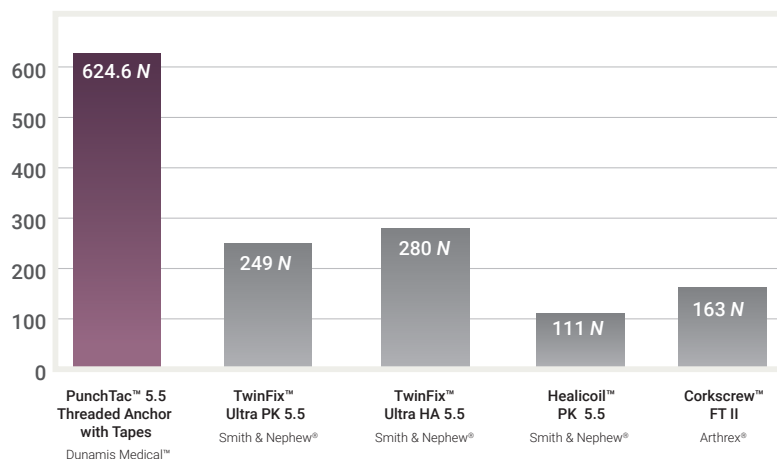

## PUNCHTac™ vPEEK Anchor System with Fiber-Connect Technology

Product No.

Box of 5

DFX0040A	5.5mm PUNCHTac™ Threaded Anchor preloaded with two 2mm Fiber-Connect™ Tape
DFX0041A	4.5mm PUNCHTac™ Threaded Anchor preloaded with two 2mm Fiber-Connect™ Tape
DFX0043A	4.75mm PUNCHTac Threaded Anchor preloaded with two 2mm Fiber-Connect™ Tape

# DUNAMIS MEDICAL™

Make the switch to better design

334.371.6280  
dunamismedical.com

<sup>1</sup> Test Data on file at Dunamis Medical.  
Some of the published competitor testing data obtained may have been performed using different density blocks.  
Vestakeep PEEK is a registered trademark of Evonik®.  
Healicoil PK, TwinFix Ultra HA and TwinFix Ultra PK are trademarks of Smith & Nephew®.  
Corkscrew FT II is a trademark of Arthrex®.  
PunchTac is a trademark of Dunamis Medical, LLC  
Listed products are covered by one or more U.S. patents and trademarks.  
Lit#: 02192019P

